



АНАЛИЗАТОР КИСЛОРОДА / АЗОТА

ELEMENTRAC ON-p

**The new ELEMENTRAC ON-p is a powerful and robust elemental analyzer for measurement of oxygen and nitrogen concentrations in inorganic materials like steel, iron, copper or ceramics. The highly sensitive NDIR and thermal conductivity detectors reliably detect element concentrations from low ppm content to high percentages.**

The innovative sample port system with pulsed chamber rinsing and vertical sample drop allow for user-friendly and comfortable analysis of rod-shaped, granular or powdery samples with a weight of up to grams.

The ELEMENTRAC ON-p elemental analyzer meets or exceeds the requirements of all relevant international standards such as ASTM E 1019 or DIN EN 3976.



[Смотреть видео](#)


**Видео о приборе**

## АНАЛИЗАТОР КИСЛОРОДА / АЗОТА ELEMENTRAC ON-P

- | Low gas consumption and high sensitivity due to closed gas system
- | Easy application of pins, powders and granules
- | Inexpensive argon as carrier gas possible
- | Short analysis time
- | Powerful impulse furnace with 8.5 kW
- | Optional Autocleaner
- | Reliable ON elemental analysis of inorganic samples like steel, non-ferrous metals, ceramics, slags, ores, etc.

**ELTRA**

## АНАЛИЗАТОР КИСЛОРОДА / АЗОТА ELEMENTRAC ON-P OPERATION AND ANALYSIS PROCESS

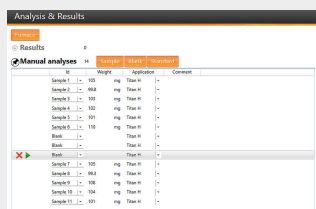

### Step 1: Logging the sample into the ELEMENTS software

The sample ID is logged into the software and the weight is automatically transferred (see step 2).



### Step 2: Weighing and introduction of sample into the port

The ELEMENTRAC ON-p analyzes volumes from a few mg up to 2 grams safely and precisely. Rod-shaped or granular samples can be applied directly. For the elemental analysis of powders, a capsule is recommended which does not have to be sealed.



### Step 3: Analysis

The empty graphite crucible is then placed on the lower electrode and the elemental analysis is started via the ELEMENTS software. The software controls all subsequent process steps.



### Step 4: Data output and export


120 to 180 seconds after the analysis has started, the measured concentrations are available for export as a report or via LIMS.

АНАЛИЗАТОР КИСЛОРОДА / АЗОТА ELEMENTRAC ON-P


## CONFIGURATIONS

The ELEMENTRAC ON -p elemental analyzer is available as a single-element analyzer for oxygen or nitrogen only, or in a multi-element configuration for measuring ON. Whereas oxygen is determined as CO<sub>2</sub> in up to two infrared cells, nitrogen is detected in its elemental form in a thermal conductivity cell.

### CUVETTE WITH VARIABLE LENGTH

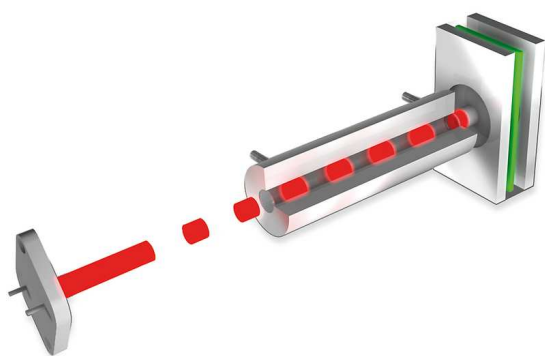

### THERMAL CONDUCTIVITY CELL WITH HIGH SENSITIVITY

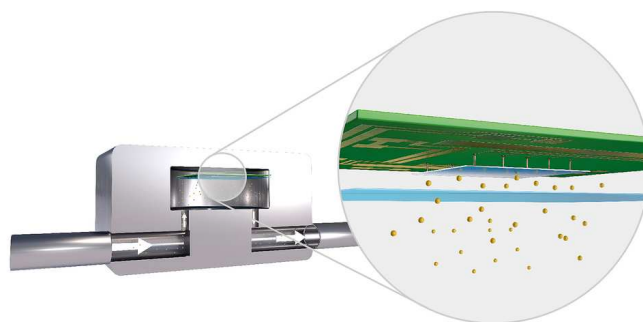

АНАЛИЗАТОР КИСЛОРОДА / АЗОТА ELEMENTRAC ON-P

## INTEGRATED STANDARD SOLUTIONS

The chemicals and filters required for operation of the elemental analyzer are arranged conveniently on the front panel and can be concealed behind a removable door during routine operation. This arrangement significantly reduces the time required for maintenance and increases user-friendliness. In addition, innovative details considerably improve the reproducibility of measurements.

### Innovative sample port & pulsed chamber flushing


The new sample port of the ON -p elemental analyzer ensures comfortable operation and reproducible measured values. Differently shaped materials like solid pieces, granules or powder in capsules can be applied up to a weight of 2000 mg, and are quickly freed from the surrounding atmosphere with the help of pulsed carrier gas flushing in the sample port. Then they drop vertically into the preheated graphite crucible for analysis.

- | Robust against dust development
- | No closing of capsules required
- | Direct application of up to 2000 mg granulate
- | Low in maintenance and wear



### Powerful catalyst

During elemental analysis in the graphite crucible, carbon monoxide (CO) is produced which is converted to carbon dioxide (CO<sub>2</sub>) in the catalyst and subsequently detected in the IR cells. The easy-to-maintain catalyst with copper oxide filling ensures complete oxidation and thus, reliable oxygen analysis even of difficult materials such as oxides.

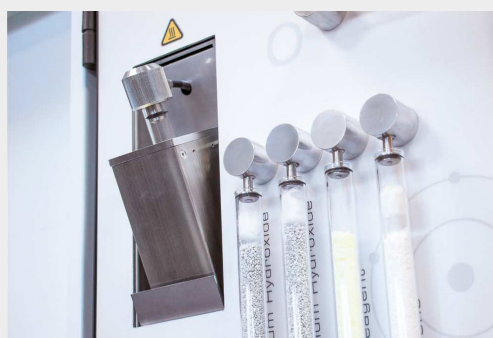

## Closed gas management

The ELEMENTRAC ONH elemental analyzer series uses a closed gas system in overpressure. This ensures that always 100% of the released sample gas is fed to the detectors which guarantees low detection limits and good reproducibility.

## АНАЛИЗАТОР КИСЛОРОДА / АЗОТА ELEMENTRAC ON-P OPTIONS

In addition to the integrated solutions of the ELEMENTRAC ON-p, further options are available to increase efficiency and extend the application range of your elemental analysis.

### Autocleaner

By melting the sample in a graphite crucible at temperatures of up to 3000 °C deposits are generated at the upper electrode and in the furnace chamber which may affect the reproducibility of ONH measurements in a negative way.

The new optional Autocleaner reliably removes these deposits, enabling precise elemental analysis even for high throughputs. Additionally, an efficient gas calibration and cleaning furnace for thorough carrier gas pre-cleaning are available for the elemental analyzer.



## ПРОГРАММНОЕ ОБЕСПЕЧЕНИЕ ELEMENTS

**The comprehensive Windows-based ELEMENTS software is an essential part of all ELEMENTRAC generation elemental analyzers.**


A central window (analysis and results) is the starting point from which all functionalities required for the daily routine are easily accessible. From here it is possible to group and export analyzed samples, or register and analyze new ones. The user may call up various subordinate functionalities like application settings, calibration, diagnosis, or status.



IOT - INTERNET OF THINGS

## THE PLATFORM FOR REMOTE ACCESS TO YOUR DEVICES

All ELTRA analyzers seamlessly integrate with the Verder Scientific IoT platform, providing enhanced functionality, seamless connectivity, and additional benefits:

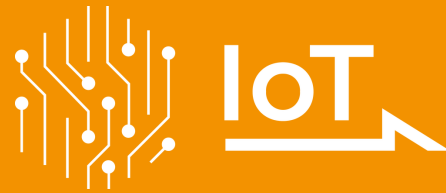

- | **Real-time Monitoring:** Gain insight into the status of your machines at any time thanks to immediate access to important data.
- | **Live Notifications:** Stay up to date on the status of your devices with instant notifications.
- | **Effortless Backup:** Whether you need to back up a single device or an entire fleet, back up your data effortlessly and minimize downtime.
- | **Automatic Software Updates:** Verder Scientific IoT keeps your device software up to date, optimizing performance and reliability.
- | **Access to Analyses Data:** ELTRA analyzers offer remote access to analysis data. This allows you to conveniently access important data while on the move.
- | **Autoloader Efficiency:** Get the most out of remote analysis preparation with our autoloader feature, which ensures uninterrupted operation and increased productivity for all instruments equipped with it.

Experience the power of the Verder Scientific IoT platform today and unlock the full potential of your ELTRA analyzers!



**FREE SOFTWARE  
DOWNLOAD**

АНАЛИЗАТОР КИСЛОРОДА / АЗОТА ELEMENTRAC ON-P

## ТИПИЧНЫЕ ОБРАЗЦЫ МАТЕРИАЛОВ

сплавы, алюминий, зола, карбиды, чугун, керамика, медь, ферросплавы, железо, металлы, руды, огнеупоры, кремний, сталь, ...



АНАЛИЗАТОР КИСЛОРОДА / АЗОТА ELEMENTRAC ON-P

## ПРИНЦИП РАБОТЫ

The measuring principle of the ELEMENTRAC ON -p elemental analyzer allows for a wide measuring range. To analyze the sample, it is weighed and placed on the sample port. Flushing with carrier gas prevents atmospheric gas (oxygen and nitrogen) from getting into the furnace.

The graphite crucible is outgassed in the impulse furnace of the analyzer to reduce possible contaminations (e.g. residual hydrogen). After a stabilization phase the sample is dropped into the crucible and melts. Carbon monoxide is produced by the reaction of carbon in the graphite crucible and oxygen of the material. Nitrogen and hydrogen are released in its elemental form. The carrier gas (helium) and sample gasses pass through a filter before entering a copper oxide catalyst which converts the CO to CO<sub>2</sub>.

The CO<sub>2</sub> is then measured by infrared cells to determine the oxygen content. CO<sub>2</sub> and water are removed chemically and the nitrogen content is measured in the thermal conductivity cell. As an option the less expensive Argon can be used in the elemental analyzer to determinate the oxygen and nitrogen content.

АНАЛИЗАТОР КИСЛОРОДА / АЗОТА ELEMENTRAC ON-P

## ХАРАКТЕРИСТИКИ

|                                       |                                                                                                                                 |
|---------------------------------------|---------------------------------------------------------------------------------------------------------------------------------|
| <b>Измеряемые элементы</b>            | азот, кислород                                                                                                                  |
| <b>Материал проб</b>                  | неорганические материалы                                                                                                        |
| <b>Направление ввода пробы в печь</b> | вертикальная                                                                                                                    |
| <b>Подача пробы в печь</b>            | графитовый тигель                                                                                                               |
| <b>Область применения</b>             | керамика, разработки / электроника, сталь / металлургия                                                                         |
| <b>Печь</b>                           | импульсная печь (максимальная мощность 8,5 кВт*) с температурой до 3.000 °C                                                     |
| <b>Принцип детектирования</b>         | метод инфракрасной абсорбции для кислорода, метод теплопроводности газов для азота                                              |
| <b>Типичное время анализа</b>         | 120 - 180 с                                                                                                                     |
| <b>Требуемые реактивы</b>             | гидроксид натрия, оксид меди, перхлорат магния                                                                                  |
| <b>Требуемый газ</b>                  | сжатый воздух, гелий качества 99.995% или лучше, аргон качества 99.995% или лучше (опционально), все газы с давлением 2 - 4 bar |
| <b>Требование к электропитанию</b>    | 3 фазы 400 В, 50/60 Гц, макс. 8,500 Вт                                                                                          |
| <b>Размеры (Ш x В x Г)</b>            | 57 x 77 x 63 см                                                                                                                 |
| <b>Вес</b>                            | ~ 161 кг                                                                                                                        |
| <b>Требуемое оборудование</b>         | ПК, монитор, весы (точность 0.0001г)                                                                                            |
| <b>Опции</b>                          | внешний чиллер, печь очистки газа-носителя, устройство газовой калибровки                                                       |
| <b>Ярлык АТМ</b>                      | * мощность программно ограничена на уровне 6.8 кВт                                                                              |

[www.eltra.com/onp2](http://www.eltra.com/onp2)

## ИНФОРМАЦИЯ ДЛЯ ЗАКАЗА

### ELEMENTRAC ON-P 2

(Пожалуйста, закажите отдельно управляющий компьютер, монитор, весы и расходные материалы (стартовый набор, ангидрон, гидроксид натрия, оксид меди II))

Measuring ranges at 1,000 mg sample weight (further measuring range combinations on request)

2)

88200-2201  ON-p 2 1xO 0.04 ppm – 0.04 % O

88200-2202  ON-p 2 2xN 0.04 ppm – 3 % N

88200-2203  ON-p 2 2xO 0.04 ppm – 1 % O

88200-2205  ON-p 2 2xO 0.04 ppm – 1 % O  
+ 2xN 0.04 ppm – 3 % N

### REQUIRED ACCESSORIES

#### УПРАВЛЯЮЩИЙ КОМПЬЮТЕР, МОНИТОР, ВЕСЫ

71015-1000 Computer with Intel Core i5-8400 Processor, 256 GB SSD; 8 GB RAM; Windows 10 operating system; keyboard; mouse

88400-0584 Monitor, TFT (23.8")

88400-0645 Весы (разрешение 0.0001 г)

### REQUIRED CONSUMABLES / CHEMICALS FOR FIRST OPERATIONS

88500-0018 ON-Starter-kit for 500 analyses  
(400 graphite crucibles, 50 outer graphite crucibles, 200 inner graphite crucibles, 50 g glass wool, 50 g quartz wool)

90200  Ангидрон (перхлорат магния), 454 г 1)

90210  Гидроксид натрия, 500 г 1)

90270  Реактив Шутца, 100 г, for OH-p and ONH-p 1)

90289  Оксид меди II, 100 г, for ON-p and ONH-p 1)

88600-0021 Copper oxide wire (for older ONH 2000 analyzer) 1)

## FURTHER OPTIONS AND CONSUMABLES

### ACCESSORIES (HARDWARE)

|            |                                                                                                     |
|------------|-----------------------------------------------------------------------------------------------------|
| 88200-2400 | ONH-p Autoloader (incl. autocleaner and vacuum cleaner)                                             |
| 88200-2401 | ONH-p Autocleaner (incl. vacuum cleaner)                                                            |
| 88400-0467 | Chiller (SMC, 5900 W)                                                                               |
| 88200-9000 | Carrier gas purification furnace, without filling (please order filling and quartz wool separately) |
| 72081      | Pressure regulator, 1 piece                                                                         |

### CRUCIBLES

|            |                                                                                                                                                                   |
|------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| 88400-0471 | Graphite crucibles, 400 pieces (recommended for autoloader operation)                                                                                             |
| 90190      |  Graphite crucibles, 400 pieces (for copper, brass and steel analysis)          |
| 90180      |  Inner graphite crucibles, 100 pieces (requires outer graphite crucible 90185) |
| 90185      |  Внешний графитовый тигель, 50 штук                                            |

### TIPS

|       |                                                                                                                                           |
|-------|-------------------------------------------------------------------------------------------------------------------------------------------|
| 31360 |  Graphite tip, 1 piece (for crucibles 90190 and 90185) |
|-------|-------------------------------------------------------------------------------------------------------------------------------------------|

### CAPSULES (REQUIRED FOR ANY KIND OF POWDER ANALYSIS)

|            |                                                                                                                                       |
|------------|---------------------------------------------------------------------------------------------------------------------------------------|
| 90257      |  Nickel capsules, 3.2 x 7 mm, 100 pieces           |
| 90256      |  Nickel capsules, 4.5 x 10 mm, 250 pieces          |
| 88400-0066 |  Nickel capsules, pressed, 12.5 x 5 mm, 100 pieces |
| 90252      |  Tin capsules, 5 x 18 mm, 100 pieces               |

#### BASKETS (REQUIRED FOR OXYGEN AND NITROGEN DETERMINATION IN REFRACTORIES)

90250  Никелевая корзинка, 1 г, 100 штук

88600-0012  Nickel baskets, high purity (low oxygen), 100 pieces, 1 g each

#### FLUXES (REQUIRED FOR SOME APPLICATIONS)

90251  Tin pellets, 454 g (for determination of hydrogen in titanium)

90800  Graphite, 50 g (improves oxygen determination)

90258  Nickel accelerator, 100 g (for analysis of high amount of refractories)

#### CHEMICALS (FILLINGS FOR GLASS AND QUARTZ TUBES)

88600-0028 Eltrasorb, 500g (black coloured sodium hydroxide)

90200  Ангидрон (перхлорат магния), 454 г 1)

90210  Гидроксид натрия, 500 г

90289  Оксид меди II, 100 г, for ON-p and ONH-p

90426-1001  Filling for carrier gas purification furnace (suitable for one filling, ONH series)

90330  Кварцевая вата, 50 г

90331  Стекловата, 454 г

90332  Стекловата, 50 г

92610  Tube of high vacuum grease, 35 g

## ELEMENTRAC - ADDITIONAL TOOLS

All ELEMENTRAC analyzers are equipped with a set of necessary tools

The following list provide part numbers for replacement of worn tools and some new tools to improve handling.

### SPATULAS AND TWEEZERS

|            |                                                                                     |                                                                                 |
|------------|-------------------------------------------------------------------------------------|---------------------------------------------------------------------------------|
| 88400-0476 |    | Micro spatula, 1 piece, XS size                                                 |
| 23110      |    | Spatula, 1 piece, M size                                                        |
| 23111      |    | Spatula, 1 piece, L size                                                        |
| 88400-0475 |    | Set with 6 spatula and 1 tweezers, for multiple weighing procedures             |
| 88400-0229 |   | Пинцет (160 мм), скругленный, 1 штука, for transporting pins and baskets        |
| 88400-0472 |  | Пинцет (145 мм), прямой, 1 штука, for removing samples out of the ONH-p furnace |
| 88400-0213 |  | Tongs for crucibles, 1 piece, for putting crucibles on the electrode tip        |

### TOOLS FOR STORAGE, TRANSPORTING AND WEIGHING

|            |                                                                                     |                                                              |
|------------|-------------------------------------------------------------------------------------|--------------------------------------------------------------|
| 88400-0477 |  | Weighing boat, 1 piece, for weighing and usage of granulates |
| 36121      |  | Quartz boat, 74x22x10 mm, 1 piece, for weighing pins         |

### TOOLS FOR CLEANING AND MAINTENANCE

|            |                                                                                     |                                                       |
|------------|-------------------------------------------------------------------------------------|-------------------------------------------------------|
| 27000-8007 |                                                                                     | O-ring set ONH-p (furnace)                            |
| 27000-8008 |                                                                                     | Maintenance kit ONH-p                                 |
| 27000-8009 |                                                                                     | O-ring set ONH-p                                      |
| 71010      |  | Brush, 16 mm, 1 piece, for cleaning balance from dust |

|            |                                                                                    |                                                                                    |
|------------|------------------------------------------------------------------------------------|------------------------------------------------------------------------------------|
| 88400-0500 |   | Telescope mirror, 1 piece, for inspection of upper electrode of ONH-p/ONH-2000     |
| 88400-0473 |   | Powder funnel (plastics), 1 piece, for easy filling of chemical tubes              |
| 88400-0489 |                                                                                    | Rubber plug 14x20x24 mm, 1 piece, for sealing small glass tubes like 88400-0006    |
| 88600-0027 |                                                                                    | Sodium hydroxide, Anhydron filter tube                                             |
| 71032      |   | Composite brush, 1 piece, for cleaning upper electrode of ONH-p furnace            |
| 71035      |   | Cleaning brush / furnace brush, 1 piece, for cleaning sample inlet of ONH furnaces |
| 71031      |   | Metal brush, 1 piece, for cleaning graphite tip and its holder                     |
| 88400-0504 |   | Cylinder brush, brass, for intensive cleaning of lower furnace                     |
| 88400-0501 |  | Micro brush, 1 piece, for cleaning of ONH series furnace outlet tube               |
| 61030      |                                                                                    | Allen key, 3 mm, 1 piece                                                           |
| 61040      |                                                                                    | Allen key, 4 mm, 1 piece                                                           |
| 61050      |                                                                                    | Allen key, 5 mm, 1 piece                                                           |

## CALIBRATION MATERIALS

**Calibration materials may show slight variations depending on the current lot.**

**To see the current certification please visit [www.ELTRA.com](http://www.ELTRA.com).**

### OXYGEN AND NITROGEN IN STEEL, PINS

|            |                                                                                     |                                                |
|------------|-------------------------------------------------------------------------------------|------------------------------------------------|
| 91100-1001 |  | Сталь, столбики 1 г, 100 штук, 25 – 40 ppm N   |
| 91100-1002 |  | Сталь, столбики 1 г, 100 штук, 30 – 70 ppm N   |
| 91100-1003 |  | Сталь, столбики 1 г, 100 штук, 150 – 250 ppm N |
| 91100-1005 |  | Сталь, столбики 1 г, 100 штук, 300 – 600 ppm N |
| 91100-1007 |                                                                                     | Сталь, столбики 1 г, 100 штук, 70 – 130 ppm N  |

|            |                                               |
|------------|-----------------------------------------------|
| 91100-1010 | Сталь, столбики 1 г, 100 штук, >1000 ppm N    |
| 91100-1011 | Сталь, столбики 1 г, 100 штук, 600-1000 ppm N |

#### HYDROGEN IN STEEL, PINS

|            |                                                                                   |                                              |
|------------|-----------------------------------------------------------------------------------|----------------------------------------------|
| 91400-1001 |  | Сталь, столбики 1 г, 100 штук, 0.5 – 1 ppm H |
|------------|-----------------------------------------------------------------------------------|----------------------------------------------|

|            |                                                                                   |                                              |
|------------|-----------------------------------------------------------------------------------|----------------------------------------------|
| 91400-1002 |  | Сталь, столбики 1 г, 100 штук, 1.5 – 4 ppm H |
|------------|-----------------------------------------------------------------------------------|----------------------------------------------|

#### СТАЛЬ, ШАРИКИ (H)

|       |                                                                                   |                                                     |
|-------|-----------------------------------------------------------------------------------|-----------------------------------------------------|
| 91110 |  | Steel, 100 balls, gold plated, 1 g each, >1.9 ppm H |
|-------|-----------------------------------------------------------------------------------|-----------------------------------------------------|

#### OXYGEN IN COPPER, PINS

|            |                                          |
|------------|------------------------------------------|
| 91000-1003 | Медь, столбики 1 г, 100 штук, ~200 ppm O |
|------------|------------------------------------------|

|            |                                                                                     |                                         |
|------------|-------------------------------------------------------------------------------------|-----------------------------------------|
| 91000-1004 |  | Медь, столбики 1 г, 100 штук, ~10 ppm O |
|------------|-------------------------------------------------------------------------------------|-----------------------------------------|

#### OXYGEN, NITROGEN AND HYDROGEN IN TITANIUM, PINS

|            |                                                                                     |                                                |
|------------|-------------------------------------------------------------------------------------|------------------------------------------------|
| 91205-1001 |  | Титан, столбики 0,1 г, 100 штук, 10 – 35 ppm H |
|------------|-------------------------------------------------------------------------------------|------------------------------------------------|

|            |                                                                                     |                                                |
|------------|-------------------------------------------------------------------------------------|------------------------------------------------|
| 91205-1002 |  | Титан, столбики 0,1 г, 100 штук, 20 – 70 ppm H |
|------------|-------------------------------------------------------------------------------------|------------------------------------------------|

|            |                                                                                     |                                                |
|------------|-------------------------------------------------------------------------------------|------------------------------------------------|
| 91205-1003 |  | Титан, столбики 0,1 г, 100 штук, 30 – 90 ppm H |
|------------|-------------------------------------------------------------------------------------|------------------------------------------------|

|            |                                                                                     |                                                 |
|------------|-------------------------------------------------------------------------------------|-------------------------------------------------|
| 91205-1004 |  | Титан, столбики 0,1 г, 100 штук, 60 – 120 ppm H |
|------------|-------------------------------------------------------------------------------------|-------------------------------------------------|

|            |                                                                                     |                                                  |
|------------|-------------------------------------------------------------------------------------|--------------------------------------------------|
| 91205-1005 |  | Титан, столбики 0,1 г, 100 штук, 150 – 250 ppm H |
|------------|-------------------------------------------------------------------------------------|--------------------------------------------------|

|            |                                                  |
|------------|--------------------------------------------------|
| 91205-1006 | Титан, столбики 0,1 г, 100 штук, 120 – 150 ppm H |
|------------|--------------------------------------------------|

#### HYDROGEN AND CARBON IN TITANIUM, PINS (250 MG)

|            |                                                |
|------------|------------------------------------------------|
| 91305-1001 | Titanium, 100 pins, 0.25 g each, < 50 ppm H    |
| 91305-1002 | Titanium, 100 pins, 0.25 g each, 50 -100 ppm H |
| 91305-1003 | Titanium, 100 pins, 0.25 g each, > 100 ppm H   |

**Please note: Every analyzer requires PC, monitor, balance and some consumables (crucibles, chemicals) which have to be ordered separately**